


# HOW TO USE ORIGIN SERVER


Software Versions: Origin Server 1.1

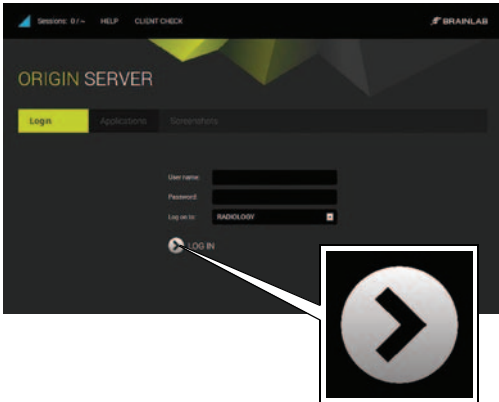

## LOG-IN


- Open your browser.
- Enter the **Origin Server** URL. Standard access follows this format:  
<https://<FQDNserver>.<hospitaldomain>>  
e.g., <https://iplanet.childrenshope.org>
- In the **Login** tab, enter your credentials and click **LOG IN**.

**NOTE:** To access additional documentation, press **Help** at the top of the **Origin Server Connection Portal**.

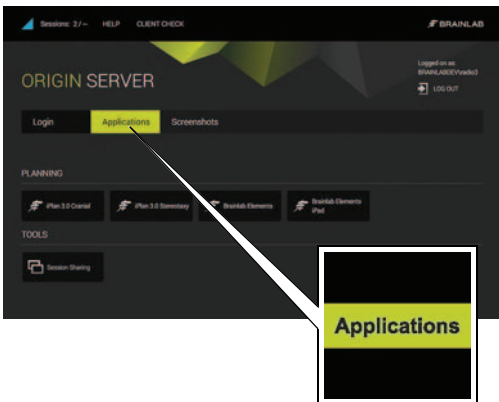

## START ORIGIN SERVER

- Start **Origin Server** by clicking on the relevant application button in the **Applications** tab.
- Perform relevant preoperative planning.

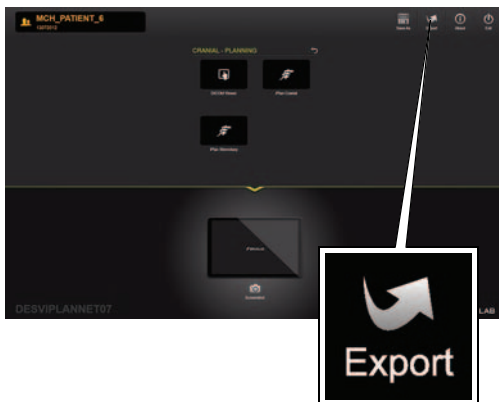

## EXPORT DATA

- Select **Export** to export your planned data.
- Wait for the data to transfer successfully.

All browsers and connection types can perform exports. However:

- USB exports must be performed using Internet Explorer in Active X mode.
- For all other browsers using RDP mode, a USB flash drive should be inserted before connection is made.
- Export using HTML5 is possible, but only when a pre-configured network folder is selected.




## EXIT ORIGIN SERVER

- Select **Done** in the application, when finished.
- Your planned data can now be used for patient treatment.
- Select **Exit**.




# ENABLE SESSION SHARING

Software Versions: Origin Server 1.1

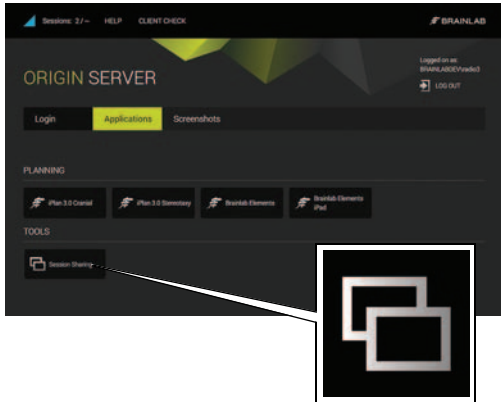

## START SESSION SHARING


- Click **Session Sharing** in the **Applications** tab. The **Session Sharing** dialog opens, listing all active sessions.

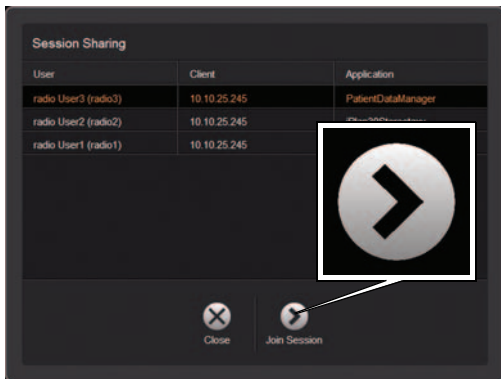

## JOIN SESSION

- Select the required session and click **Join Session**.
- Accept the session as prompted.

**NOTE:** Both the session owner and the guest user can now use the software's full functionality and see all steps performed by the other user.

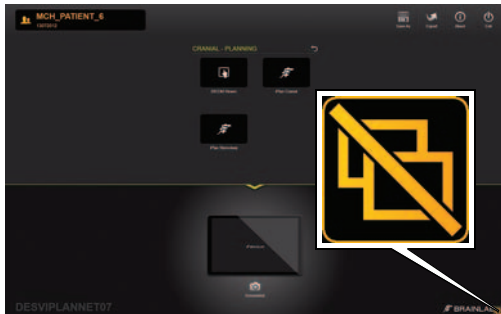

## DISCONNECT THE SESSION


- To leave the shared session, click **Disconnect** at the bottom right of the screen.

**NOTE:** When you disconnect, the application continues to run and does not close.



# ORIGIN SERVER WITH TABLETS

Software Versions: Origin Server 1.1

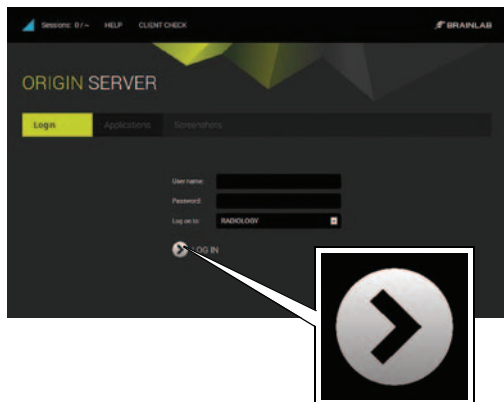


## LOG-IN

- Open your HTML5 compatible browser.
- Navigate to **Origin Server** using the web browser on the tablet device. Standard access follows this format:  
<https://<FQDNserver>.<hospitaldomain>>  
e.g., <https://iplanet.childrenshope.org>
- In the **Login** tab, enter your credentials and click **LOG IN**.

**NOTE:** Wireless access to the server is required.

**NOTE:** Brainlab recommends using Safari.

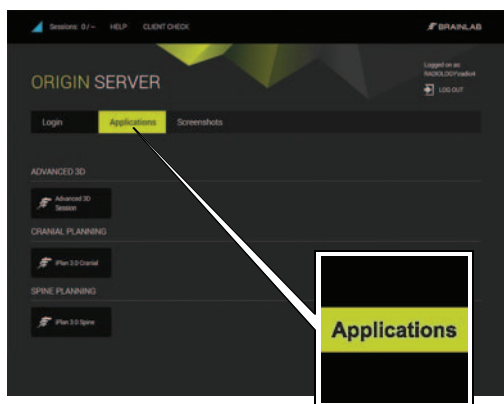

## SELECT APPLICATION

- Start e.g., **Brainlab Elements iPad**.
- Start an application by clicking on the relevant application in the **Applications** tab.
- Perform relevant preoperative planning.


**NOTE:** Do not use tablets for diagnostic purposes. Use tablets for information only.



## ADD TO HOME SCREEN - STEP 1

- Select the menu icon ①.  
Menu opens.
- Select Add to Home Screen icon ②.

**NOTE:** For optimized viewing it is recommended to view **Origin Server** in full screen mode. Adding the server URL to the Home Screen ensures that **Origin Server** opens each time in full screen mode.

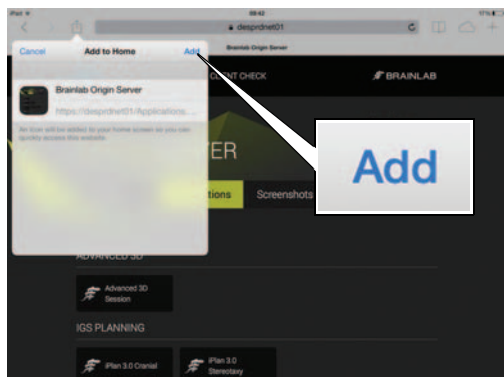


## ADD TO HOME SCREEN - STEP 2

- Select **Add** from the dialog.  
An icon is added to your home screen.



## START ORIGIN SERVER

- Start **Origin Server** by selecting the icon you created from the Home Screen.

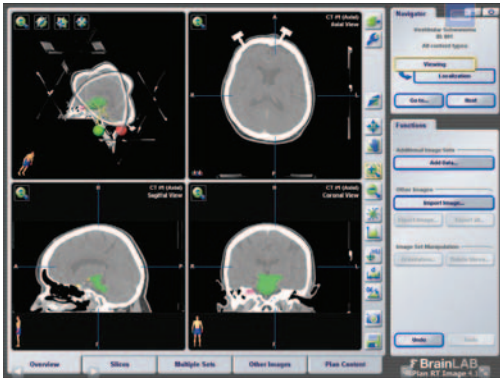


## VIEWING

- It is recommended to use full screen capable browsers.

**NOTE:** Exact viewing functionality is browser-dependent and the chosen application may not always fit exactly to the size of the screen.

**NOTE:** Do not navigate back and forth during surgery planning.

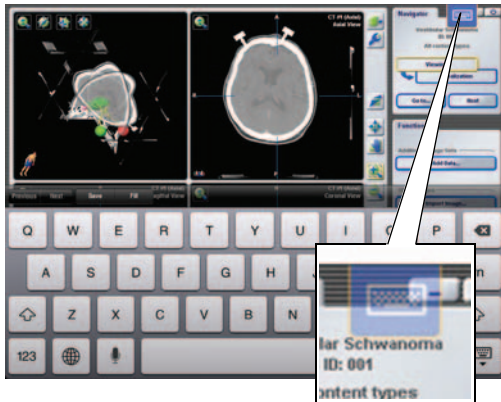

## ACTIVATE THE KEYBOARD

- Activate the keyboard by selecting the keyboard icon from the overlay menu.

**NOTE:** Exact functionality is browser-dependent.

### COPYRIGHT:

This guide contains proprietary information protected by copyright. No part of this guide may be reproduced or translated without the express written permission of Brainlab.

Document Revision: 1.0  
Article Number: 60917-10EN

### LIABILITY:

This guide is subject to change without notice and does not represent a commitment on the part of Brainlab.

For further information, please refer to the "Limitations of Liability" section in the Brainlab Standard Terms and Conditions of Sale.

